

TRANSPORT SPECIMENS SAFELY



Zorb Sheets | Blood & Body Fluid Absorbents


For diagnostic specimens, infectious substances, and dangerous goods

- Available in multiple sizes
- Fast-acting and easy-to-use
- Lightweight (space saving)
- Prevents leaking (transforms fluid into semi-solid gel-like mass)

How Zorb Sheets Work

The tissue-based material allows fluid to immediately activate the super absorbent polymers to solidify and transform the contents into a semi-solid (gel) mass. In the event of a spill or leak, it creates an extra protection barrier for workers' safety. (These high quality granules have the unique ability to hold large amounts of absorbed water-based fluid, even against pressure.)

Item No.: 44002
Size: 12" x 12"
Case Qty: 500
Max. Blood/Urine/Saline Absorbancy: Up to 275cc
Max. Water Absorbancy: Up to 720cc

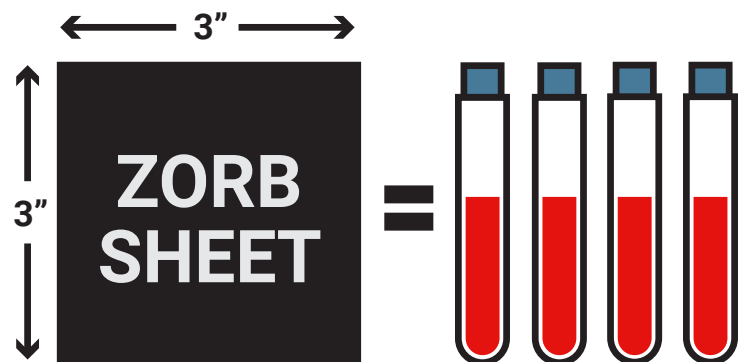

Just one 3"x3" Zorb Sheet absorbs 20ml of blood; that's more than four 5ml samples.

 **PharmaSystems®**

PharmaSystems.com
151 Telsons Road
Markham Ontario
Canada L3R 1E7