

SE46HCRxGV

Pharmaceutical /
Laboratory Refrigerator



Embedded Digital
Data Logger



Specifications

Width*	Depth**	Height***	Shelves	Interior Cubic Volume	Ship Weight	Ship Cube	Cooling Capacity	Electrical
47 1/2" 1207 mm	31" 788 mm	78" 1982 mm	10	37.21 ft ³	415 lbs 188 kg	78 cu. ft 2.21 m ³	3/4 HP	115V/7.42 Amp/60 Hz - Single Phase - NEMA 5-15P, Dedicated Circuit

** Add 1" (26mm) for Power Cord. *** Add 3/8" (10mm) for Threaded-In Levelers. Millimeters rounded up.

Premium Features



Automatic Stainless Steel
Condensate Control System



Cassette®
Refrigeration



Free & Clear
Condenser Coil



Glycol Bottle
Probe



Embedded DDL



ProductTemp™



Variable speed



Digitally Controlled



Front Breathing
Refrigeration




Priority Air-Flow

SE46HCRxGV: Designed to enhance pharmaceutical storage and display experience to the highest standards of quality and performance.

Refrigeration:

- Self contained Cassette® forced-air refrigeration employs R290 hydrocarbon Natural Refrigerant within a sealed capillary tube system. Cassette® refrigeration is bottom-mounted and removable from the front for ease of maintenance. Refrigeration may be exchanged or upgraded without mechanic. Certified to NSF-456 for temperature holding (35.6°F to 46.4°F / 2°C to 8°C) and recovery.
- Equipped with Energy Efficient Digital Controller featuring visual/audible product protection alarm.
- Create own heading named "Monitorings and Alerts".
- Priority, directed air flow provides optimum cooling and energy efficiency within cabinet and exclusively at front of unit for optimized heat exchange ventilation.
- Finless, Free&Clear™ Condenser Coil maximizes air-flow over the compressor, reduces energy consumption and is non-clogging.
- Automatic, stainless steel condensate control system.
- Technologically advanced ProductTemp™ provides exterior product temperature display. Displays to a resolution of 0.1°.
- High quality EC motors providing streamlined airflow, quieter and yet robust refrigeration.
- Variable speed compressor for extended operational life, improved energy efficiency and reduced pull down time on demand to keep tight temperature fluctuations.

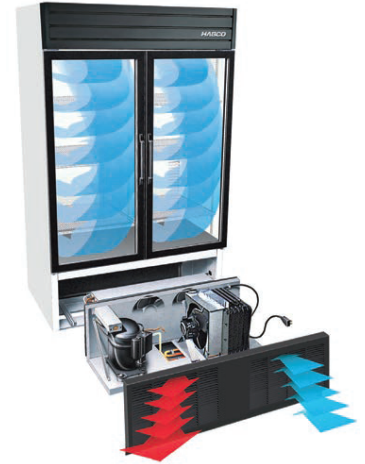

Cabinet:

- Exterior – durable appliance white.
- Interior – NSF approved white aluminum liner with stainless steel floor.
- Equipped with a Probe Port, which allows for the insertion of a Third-Party Temperature Probe.
- Ergonomic raised interior floor.
- Space Saving Design – air ventilation exclusively from the front of unit.
- Insulation – cabinet is foamed-in-place with eco-conscious, polyurethane insulation.
- Heavy duty thread-in leveler legs.

Doors:

- Self Closing Swing Doors feature hold-open position, spring back protection and are reversible on-location.
- Door Lock is standard.
- Magnetic gaskets are single piece bead-mount construction easily removable for cleaning.
- Heat Reflecting HABCO Insulated Safety Glass with black extruded frame.
- Integrated, injection molded pocket handle(s).

Lighting:



- Energy Efficient LED Lighting System for sparkling Illuminated Interior Display, on/off Switch.

Storage:

- Ten (10) PU coated white shelves 21 1/4"W x 20 1/2"D (540mm x 521mm), are NSF approved.
- Heavy duty aluminum shelving pilasters will not peel or chip adjustment of 1/2" (13mm) increments.

Electrical:

- Factory power cord installed for connection to a sole-use (dedicated) 115V/60Hz/1-Phase, 15 Amp circuit.

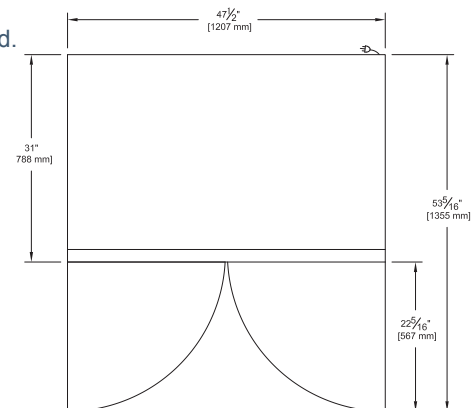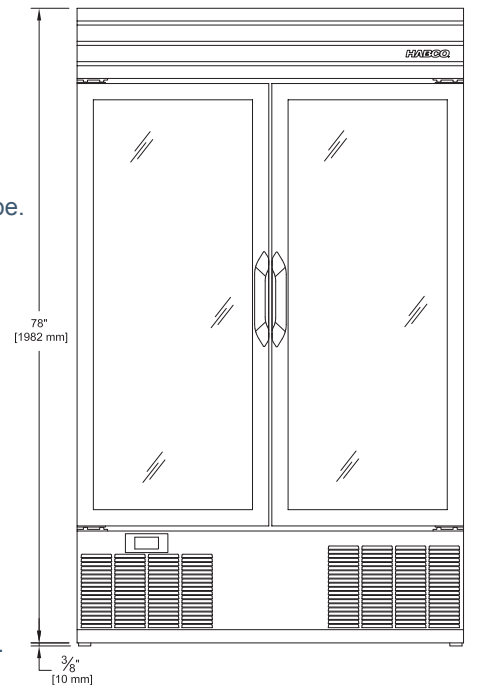

HABCO MANUFACTURING INC. 501 GORDON BAKER ROAD, TORONTO, ONTARIO, CANADA, M2H 2S6
Toll Free CAN:1-800-387-5398 USA:1-800-448-0244 Fax:416-491-6982 www.habcomfg.com